

# Agilent xCELLigence Immunotherapy IMT Assay (anti-CD40) Tethering Kit and Sample Kit

Recommended Use: Monitoring liquid tumor cell killing

For Research Use Only. Not for use in diagnostic procedures.

## Kit components

**Table 1.** List of kit components.

Product	Part Number
IMT Assay (anti-CD40) Tethering Kit (up to 6 Plates)	8100005
125x Tethering Reagent (anti-CD40) 250 $\mu$ L	
10x Tethering Buffer, 10 mL	
Cytolysis Reagent, 1.5 mL	
IMT Assay (anti-CD40) Sample Kit (up to 2 Plates)	8100006
125x Tethering Reagent (anti-CD40) 90 $\mu$ L	
10x Tethering Buffer, 10 mL	
Cytolysis Reagent, 1.5 mL	

## Storage conditions

Tethering Reagent and Cytolysis Reagent should be stored undiluted at 4 °C.

**Do not freeze.**

The 10x Tethering Buffer and E-Plates can be stored at 25 °C.

## Additional materials and reagents required

The following are required but are not supplied with these kits:

- Tissue culture media
- Tissue culture flasks/plates
- PBS without calcium and magnesium
- Sterile tissue culture grade water
- Agilent xCELLigence instrument with RTCA Software Pro
- Liquid tumor target cells
- Immune effector cells
- Small molecule drugs and/or biological effector molecules (antibodies, etc.)
- E-Plate


## Introduction

Agilent xCELLigence real-time cell analysis (RTCA) instruments use gold biosensors embedded in the bottom of microtiter plate wells to noninvasively monitor the status of adherent cells using the principle of cellular impedance. Adherent cells act as insulators by impeding the flow of an alternating microampere electric current between the biosensors. This impedance signal is measured automatically, at a frequency defined by the user (every 10 seconds, once per hour, etc.), and provides a sensitive readout of cell number, cell size/shape, and cell-substrate attachment strength.

Over the past decade xCELLigence has been used to study immune cell-mediated killing of adherent cancer cells. This IMT Assay Kit expands the repertoire of xCELLigence to include monitoring immune cell-mediated killing of nonadherent liquid cancers. Below, the xCELLigence RTCA assay principle and workflow are first described using adherent target cells as an example. Following a general overview of how the xCELLigence RTCA assay has been adapted for monitoring liquid tumor killing, a detailed protocol is provided.

### xCELLigence IMT assay principle

In contrast to adherent cancer cell targets, immune effector cells are nonadherent and therefore produce a minimal impedance signal. Because of this, when adherent cancer cells are treated with various effectors (NK cells, T cells, CAR-T cells, oncolytic viruses, checkpoint inhibitors, bispecific antibodies, BiTEs, etc.) it is possible to selectively monitor the kinetics of cancer cell destruction in real time. The major distinguishing features of this technology include enhanced sensitivity, the preclusion of labels, simple workflow and continuous kinetic measurement of cancer cell health/behavior. The xCELLigence immunotherapy workflow is outlined in Figure 1.

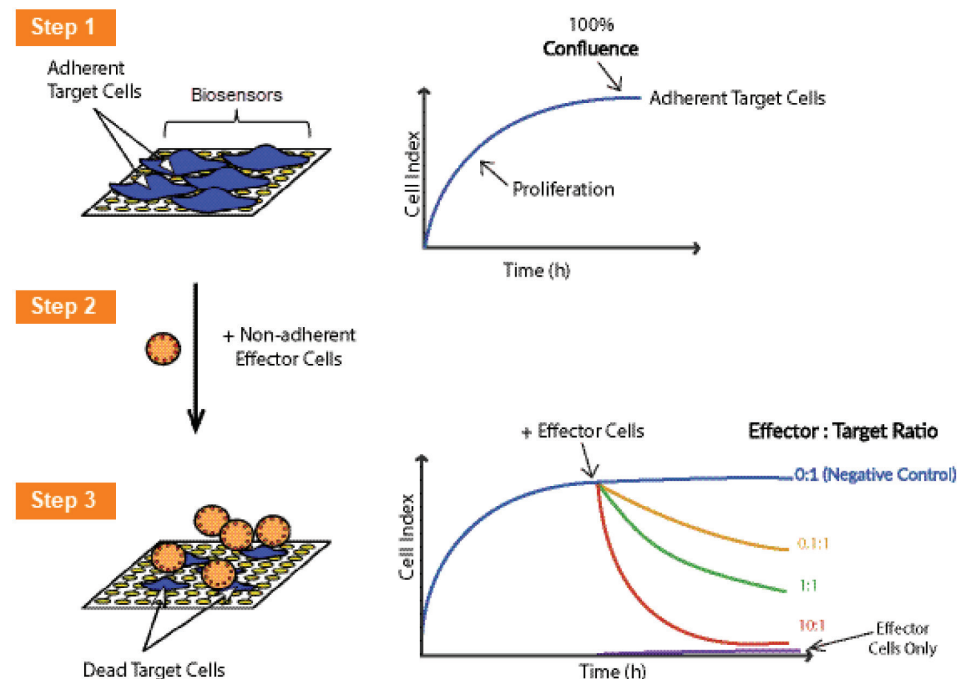

**Figure 1.** Monitoring immune cell-mediated killing of adherent cancer cells in real time using Agilent xCELLigence.


**Step 1:** Adherent target cells (i.e. tumor cells) are first seeded in the wells of an electronic microtiter plate (E-Plate). Adhesion of cells to the gold biosensors impedes the flow of electric current between them. This impedance value, plotted as a unitless parameter called "Cell Index", increases as cells proliferate and then plateaus as cells approach 100% confluence.

**Step 2:** When added subsequently, nonadherent effector cells (i.e. immune cells) in suspension do not cause changes in impedance (due to lack of adherence to the gold biosensors).

**Step 3:** If effector cells induce the destruction of the target adherent tumor cells, the corresponding cytolytic activity can be sensitively and precisely detected. The continuous acquisition of impedance data for each well of an E-Plate enables the generation of real-time killing curves for multiple conditions simultaneously.

### **Using xCELLigence to monitor liquid tumor cell killing**

The dozens of peer-reviewed studies published over the past decade have established the utility of xCELLigence RTCA for probing the efficacy of immunotherapies targeting solid/adherent cancers. However, ~10% of all cancers are liquid in nature, and therefore nonadherent, and cannot be directly monitored by the standard impedance assay. Liquid cancers (B cell cancers in particular) are prominent immunotherapy targets at present because they are readily accessible within the bloodstream and are not confounded by the microenvironment complexities/heterogeneities associated with solid tumors. To help accelerate research in this area, Agilent developed this xCELLigence IMT Assay Kit focused on liquid tumor cell killing. For this purpose, the wells of Agilent E-plates are precoated with an anti-CD40 Tethering Reagent, enabling CD40+ liquid tumor cells to be immobilized on the plate bottom (Figure 2A). Whereas Tethering Reagent immobilized liquid tumor cells generate a robust impedance signal and proliferate to the point of confluence (resulting in a plateaued impedance signal), the growth of untethered cells is undetectable (Figure 2B). Importantly, with or without anti-CD40 coating of the wells, effector cells such as the NK-92 cells used here produce minimal signal on their own (Figure 2B). Addition of NK-92 cells on top of immobilized liquid tumor cells results in target cell death in a dose-dependent manner (Figure 2C). The tethering and killing behaviors seen in Figure 2B and Figure 2C have been observed in all three of the B cell lymphoma cell lines tested (Daudi, Raji, and Ramos), for multiple effector cell types (NK, T, CAR-T), and for combination therapies (CAR-T and checkpoint inhibitors, etc.).

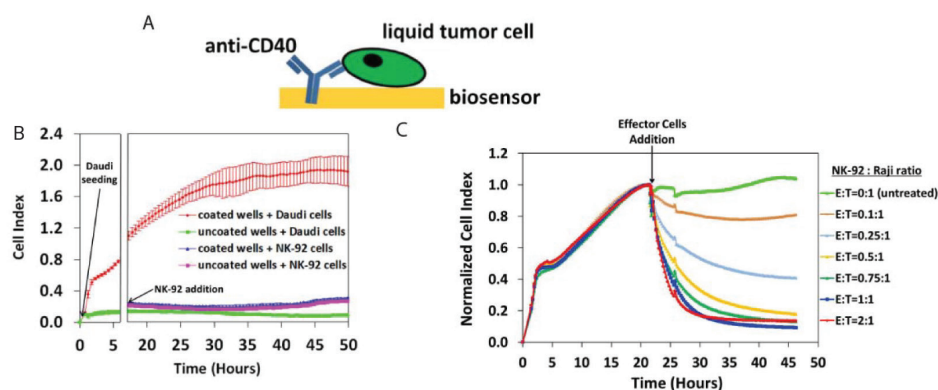

**Figure 2.** The Agilent xCELLigence IMT Assay Kit to monitor liquid tumor cell killing. (A) Precoating the wells of Agilent E-Plates with anti-CD40 enables liquid tumor cells to proliferate on and be detected by the biosensors. (B) Controls showing the selective proliferation of Daudi B lymphoblasts on biosensors coated with anti-CD40. As expected, with or without anti-CD40 coating nonadherent NK-92 effector cells produce minimal signal. Error bars are standard deviation. (C) Raji target cells were seeded at 60,000 cells/well. The efficiency with which target cells are killed is dependent on the number of effector cells added per well.

## Precautions

See the [Safety Data Sheet](#)

## Protocol

### Overview

The below protocol has specifically been optimized for using the xCELLigence RTCA SP or MP instruments to monitor effector cell mediated killing of B lymphoblasts (Daudi, Raji, and Ramos cell lines were used in the optimization of this protocol) which have been immobilized on the well bottoms of an E-Plate View 96 using the Tethering Reagent (anti-human CD40 antibody) and Tethering Buffer provided with this kit. Assay conditions may require additional optimization if different cell lines or antibodies are used.

### Reagents and equipment

The following reagents and equipment were used in developing/optimizing this protocol:

#### Reagents

- **Cells:** Daudi Burkitt's lymphoma B lymphoblasts (part number CCL-213), Raji Burkitt's lymphoma B lymphoblasts (part number CCL-86), and NK-92 non-Hodgkin's lymphoma cytotoxic NK cells (part number CRL-2407) were purchased from ATCC.
- **Cell Culture Medium:** RPMI 1640 + 10% FBS + 1% penicillin/streptomycin
- **Interleukin 2 (IL-2):** Peprotech (part number AF-200-02). Working concentration is 200 units/mL for NK-92 cells.

## Equipment

- **E-Plate View 96:** The media volumes and cell numbers used in the protocol below are compatible with usage of the lower throughput E-Plate 16 because the well size in the Agilent E-Plate View 96 is identical to the well size in the E-Plate 16. the Agilent E-Plate 384 is also compatible with this IMT assay, but cell numbers/media volumes would need to be optimized for the different well dimensions in these plates.
- **xCELLigence RTCA SP or MP instruments:** While all of the xCELLigence RTCA instruments can be used for this IMT assay, the cell numbers and media volumes listed in this protocol are optimized for the dimensions of the wells in the E-Plate View 96, which is only compatible with the xCELLigence RTCA SP and MP instruments.

## Workflow summary

Allowing an initial day for target cells to attach to and proliferate in the E-Plate View 96 wells, this assay protocol has been developed for continuous monitoring of cell killing over the course of 4 subsequent days.

### Day 1

1. Coat an E-Plate View 96 with diluted Tethering Reagent for 3 hours at room temperature. Alternatively, coat the E-Plate View 96 with diluted Tethering Reagent at 4 °C overnight. Wash, add medium, and take a background measurement.
2. Prepare target cells and add to the coated wells of the E-Plate View 96.
3. Incubate the E-Plate View 96 at room temperature for at least 30 minutes to let the cells settle.
4. Load the E-Plate View 96 into the xCELLigence RTCA instrument and start data acquisition to monitor target cell attachment and proliferation.

### Day 2

5. Prepare the effector cells.
6. Pause xCELLigence data acquisition; remove the E-Plate View 96 from the instrument and place inside a hood; remove media from the wells; add effector cells to the E-Plate View 96.
7. Load the E-Plate View 96 back into the xCELLigence RTCA instrument and start data acquisition to monitor effector cell-mediated killing of immobilized target cells.

### Days 2 to 5

8. Continue data acquisition for as long as desired.
9. Analyze the data.

## Detailed protocol

### Coating the E-Plate View 96 with Tethering Reagent (3 hours)

1. Tethering Buffer is provided as 10 mL of a 10x solution. Dilute to 1x by adding 90 mL of sterile tissue culture grade water to the bottle.
2. The Tethering Reagent (anti-CD40) is provided at a concentration of 500 µg/mL. Using 1x Tethering Buffer, dilute the Tethering Reagent 125-fold to a final concentration of 4 µg/mL.
3. Coat the E-Plate View 96 by adding 50 µL of the 4 µg/mL Tethering Reagent to each well.
4. Incubate at room temperature for 3 hours. Alternatively, incubate the E-Plate View 96 at 4 °C overnight.

**CAUTION:** The extent of Tethering Reagent adhesion to well bottoms may vary as a function of incubation temperature and duration. To ensure reproducibility between experiments be sure to maintain consistency in the parameters of this coating step.

### Washing and background measurement (5 minutes)

5. Remove Tethering Reagent.
6. Gently wash wells twice with 200 µL of PBS.
7. To each well, add 50 µL of complete medium (or the medium to be used during the assay).
8. Place the E-Plate View 96 back into the instrument and let it equilibrate at 37 °C for 10 minutes before taking a background measurement.

### Target cell seeding (45 minutes)

9. Determine target cell concentration and seed 40,000 to 60,000 cells in a volume of 100 µL/well. With the medium present in the wells (from the background measurement step), this gives a final volume of 150 µL/well.

**CAUTION:** The number of cells used in this IMT assay will ultimately depend on the cell type being used. Conducting preliminary experiments to determine the optimal cell number is highly recommended.

10. Incubate the E-Plate View 96 at room temperature for 30 minutes to facilitate uniform immobilization of target cells on plate bottom.

**CAUTION:** Failure to perform this step can result in large well to well variation in the cell seeding density/pattern (which can affect the impedance signal) for the following reason: Immediate warming to 37 °C can cause convection currents to form within the well, pushing cells to the well perimeter and resulting in an uneven distribution of cells on the impedance biosensors.

This well-known phenomenon is not unique to Agilent E-Plates; it occurs in all microtiter plates.

### **Monitoring target cell adhesion and proliferation (24 hours)**

11. Place the E-Plate View 96 back into the xCELLigence instrument and initiate data acquisition. The purpose here is simply to monitor target cell adhesion and proliferation prior to initiating cell killing. Accordingly, measuring impedance every 15 minutes for 24 hours is recommended.

### **Effector cell addition (15 minutes)**

12. Effector cells (NK-92 as an example here) are typically grown in the presence of interleukin 2 (IL-2) at a concentration of 200 units/mL (NK-92 cells). To avoid IL-2-mediated stimulation of target cells during the IMT assay, pellet effector cells and resuspend in media without IL-2. Determine effector cell concentration and prepare serial dilutions.
13. Press the skip button to fast-forward through the remaining sweeps in the currently running step of data acquisition; do not use the pause button. Remove E-Plate View 96 from instrument and place it in the hood.
14. Gently aspirate 50  $\mu$ L of medium from each well, leaving 100  $\mu$ L in each well.
15. Add 100  $\mu$ L of the NK-92 dilutions per well, being careful to not disturb or remove the immobilized target cells. Effector:Target cell ratios will be defined by the density of the effector cells in the serial dilutions.
16. Incubate the E-Plate View 96 at room temperature for 30 minutes to facilitate uniform distribution of effector cells on top of immobilized target cells.

### **Cytolysis positive control (10 minutes)**

17. For a positive control demonstrating rapid and complete target cell killing, perform the same steps used for effector cell addition, but in place of effector cells, add 100  $\mu$ L of Cytolysis Reagent that has been diluted 10-fold using culture medium. This detergent based reagent will cause the Cell Index to drop to  $\sim$ 0 within minutes.

### **Monitoring the killing of immobilized target cells (assay duration defined by the user)**

18. Place the E-Plate View 96 back into xCELLigence instrument and initiate data acquisition. The purpose here is to monitor the killing of the target cells immobilized on the bottom of the E-Plate View 96 wells. Impedance measurements every 15 minutes are recommended. The total duration of data acquisition should be programed to be longer than deemed necessary; data acquisition can be terminated whenever appropriate.

### **Data analysis**

The full functionality of the RTCA data analysis and representation software is reviewed in detail in the RTCA Software Pro manual.

## Troubleshooting

Potential problems and corresponding solutions are listed below:

### **There is large variation between replicate wells**

- It is important to perform the full 30-minute room temperature incubation after seeding both target cells and effector cells. This ensures even distribution of the cells across the biosensors.
- Ensure cells are well suspended prior to seeding into the E-Plate View 96 for consistent cell seeding densities between wells.
- The reverse pipetting technique is effective for preventing bubble formation and assuring consistent volume dispensing.

### **No signal is detected in wells expected to have a signal from target cell immobilization**

- Verify that target cells have been immobilized on the E-Plate View 96 well bottoms by replacing cultured medium (which may contain uncaptured cells) with fresh medium and then examining wells under a microscope.
- Ensure the use of the Tethering Reagent (anti-CD40) and Tethering Buffer provided with this kit at the recommended dilutions. Reagents have been optimized for this kit.
- Verify CD40 expression on the surface of the target cells using flow cytometry.

[www.agilent.com/chem/xCELLigence](http://www.agilent.com/chem/xCELLigence)

RA44349.2464699074

This information is subject to change without notice.

© Agilent Technologies, Inc. 2020, 2021  
Printed in the USA, July 14, 2021  
5994-3238EN